

HAIRMAX®
The science of hair growth.




HairMax LaserComb®

Prima 7

QUICK START GUIDE


The HairMax LaserComb® is a very effective device when used correctly. It is manufactured using precision components and should last 10-15 years under normal usage. Please take time to read the entire User Manual before operating the LaserComb.

Prima 7 LaserComb Model

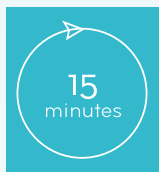

Use the Prima 7 for 15 minutes, 3 times a week. The treatments should be spaced every other day, for example Monday, Wednesday and Friday.


It is recommended, but not essential, that the LaserComb be used on a clean scalp after washing. Your hair can be wet, towel dried or dry. The LaserComb should be used prior to applying any topical treatments, gels, sprays or concealing products.

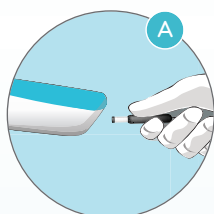

- **Charge your LaserComb** by placing the power cord directly into the handle end of the LaserComb (Figure A).

For initial charging, allow 2-3 hours to fully charge your LaserComb.


The indicator light on the power button will flash red during the charging process. Once charging is complete, the light will stop flashing and remain on.

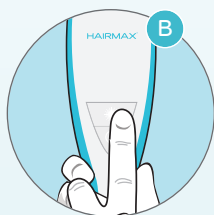

- **Turn on the LaserComb** by pressing the power button and hold until the power button light turns on.

Next, press the laser button until it lights up. That means the LaserComb is ready for use (Figure B).

- **Before beginning treatment**, be sure to test the LaserComb for sensitivity (see page 5 of your User Manual). To start your treatment place the LaserComb flat on your scalp so that both rows of teeth are touching the scalp (Figure C).

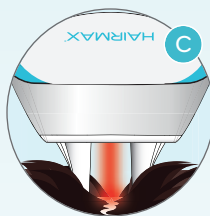

- **Start at your hairline**. Move the LaserComb from front to back of scalp in a slow combing motion, stopping every ½ inch (1cm) and pausing for 4 seconds (Figure D). Do not lift the LaserComb off your scalp while combing from front to back.



- **Lift the LaserComb off your scalp** and start at a new position at the hair line. Continue in this manner from front to back and then side to side (Figure E) until you cover your entire scalp with multiple passes. When the LaserComb is on it will beep and vibrate every four (4) seconds. This means that it is time to move the LaserComb to the next area of the scalp. Repeat these passes for the total treatment of 15 minutes.



- **If you have long hair**, use your fingers to lift your hair out of the way while using the LaserComb to be sure that the laser energy reaches your scalp (Figure F).



- **After treatment**, turn off the LaserComb by pressing and holding the power button (lower button) until the indicator light turns off (Figure G).



Please Note:

- The total time required for treatment with the HairMax Prima 7 LaserComb is 15 minutes.
- It is important to move the LaserComb very slowly. Remember, you are bathing the scalp with laser light. If you move it too fast, the treatment will not achieve the optimal results to promote hair growth.
- While there are no negative effects from overuse, it is recommended that you use the LaserComb according to directions.
- The LaserComb has advanced electronics to make sure that the device is working properly. If for some reason you are experiencing difficulties with your device, before sending it back to us for repair, please refer to the Trouble Shooting section of your User Manual.

For further assistance please call +1-561-417-0200.



REGISTER

To validate your warranty please register online at www.hairmax.com/warranty or call Customer Service toll-free at 1-866-527-3726

CONTACT US AT LEXINGTON INTERNATIONAL, LLC:

Corporate Office - Customer Service - Warranty Repairs
Lexington Intl., LLC
1040 Holland Drive
Boca Raton, Florida 33487 USA

United States 1-866-527-3726
1-561-417-0200
United Kingdom 0808-134-9868
Australia 1-800-937-358

Patent www.HairMax.com/patents

© Copyright 2017 Lexington Intl., LLC.
All rights reserved.