


IS HAIRMAX® RIGHT FOR YOU?

HairMax is most effective in women and men in early to moderate stages of hair loss and thinning. If your hair loss looks like the images below, HairMax is recommended for you.

FEMALE



MALE

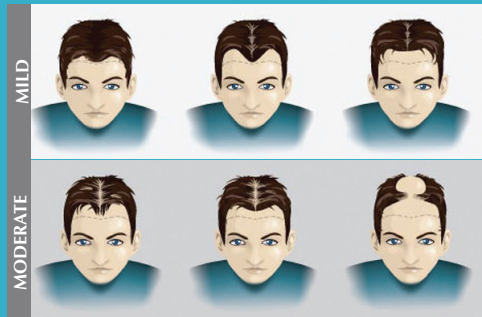

DOCTOR RECOMMENDED

Physicians worldwide recommend HairMax Laser Light Therapy to their patients. HairMax is supported by a Medical Advisory Board of established leaders in the fields of dermatology and hair restoration.

FREQUENTLY ASKED QUESTIONS

How often do I need to use my LaserBand 41?

It is recommended to use your LaserBand 41 three times a week for as little as 3 minutes, with at least a day between each treatment.

Will results improve if I use my LaserBand 41 more than 3 times a week?

Three times a week has been found to be the optimal treatment time. Using your LaserBand 41 more than 3 times a week has not shown any additional benefits.

Should I use my LaserBand 41 on wet or dry hair?

You can use your device on dry or damp hair. It is important that your scalp be clean and free of styling products. A build-up of topicals such as oils, gels, sprays, etc. can block the laser light from effectively reaching the hair follicles.

Do I have to use any special shampoos or products with my LaserBand 41?

No. However, we recommend using the HairMax® for den•sity collection. These paraben-free, sulfate-free, low pH products are gentle on your thinning hair and complement the laser treatment.

HAIRMAX THINNING HAIR CARE

From roots to ends HairMax brings new life to your thinning hair with the den•sity collection, a full line of specially formulated hair care products including shampoo, conditioner, scalp treatments, supplements and hair fibers.




Lexington International, LLC
1040 Holland Dr, Boca Raton, FL 33487
Tel. +1.561.417.0200 | U.S. Toll Free 1.800.9REGROW
FAX. +1.561.892.0747

www.HairMax.com

LIT-CONBROLEX-LB41 00 12/16 © Copyright 2017 Lexington Intl., LLC. All rights reserved.



NEXT GENERATION HAIR GROWTH LASER DEVICE

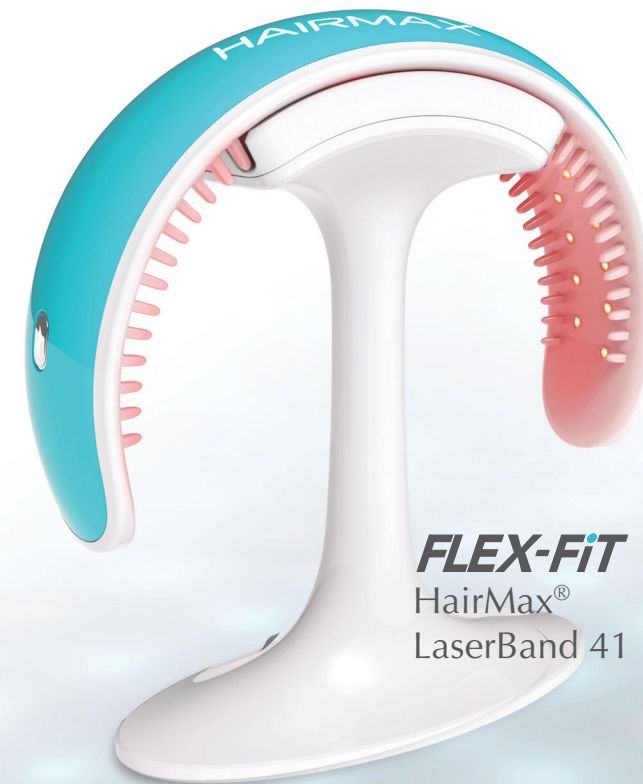


FLEX-FIT
HairMax®
LaserBand 41

HAIRMAX®
The science of hair growth.™

HOW LASER LIGHT THERAPY PROMOTES HAIR GROWTH

HairMax® works through a process called Photobiostimulation, delivering therapeutic, nourishing light energy to your hair follicles. This boost of laser light revitalizes and enlivens weakened hair follicles, restores your hair's natural growth cycle, reverses the thinning process and stimulates hair growth with increased density, fullness and vibrancy.

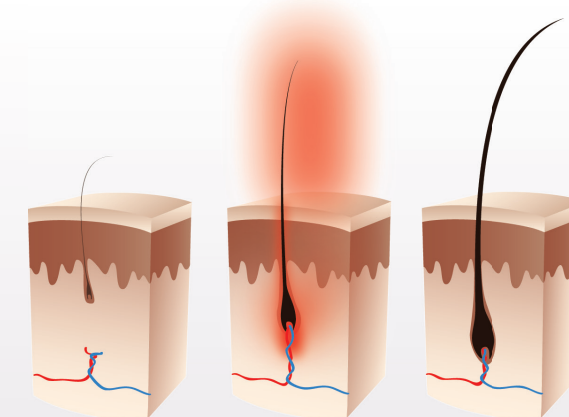


LASER PHOTO THERAPY FOR HAIR GROWTH

CLINICALLY PROVEN

HairMax laser devices have been proven safe and effective in 7 clinical studies with a total of 460 men and women. Participants experienced an average increase in hair count of 129 additional new hairs per sq. in.

90% participants in clinical studies experienced significant hair growth

Next Gen Laser Light Hair Growth Device

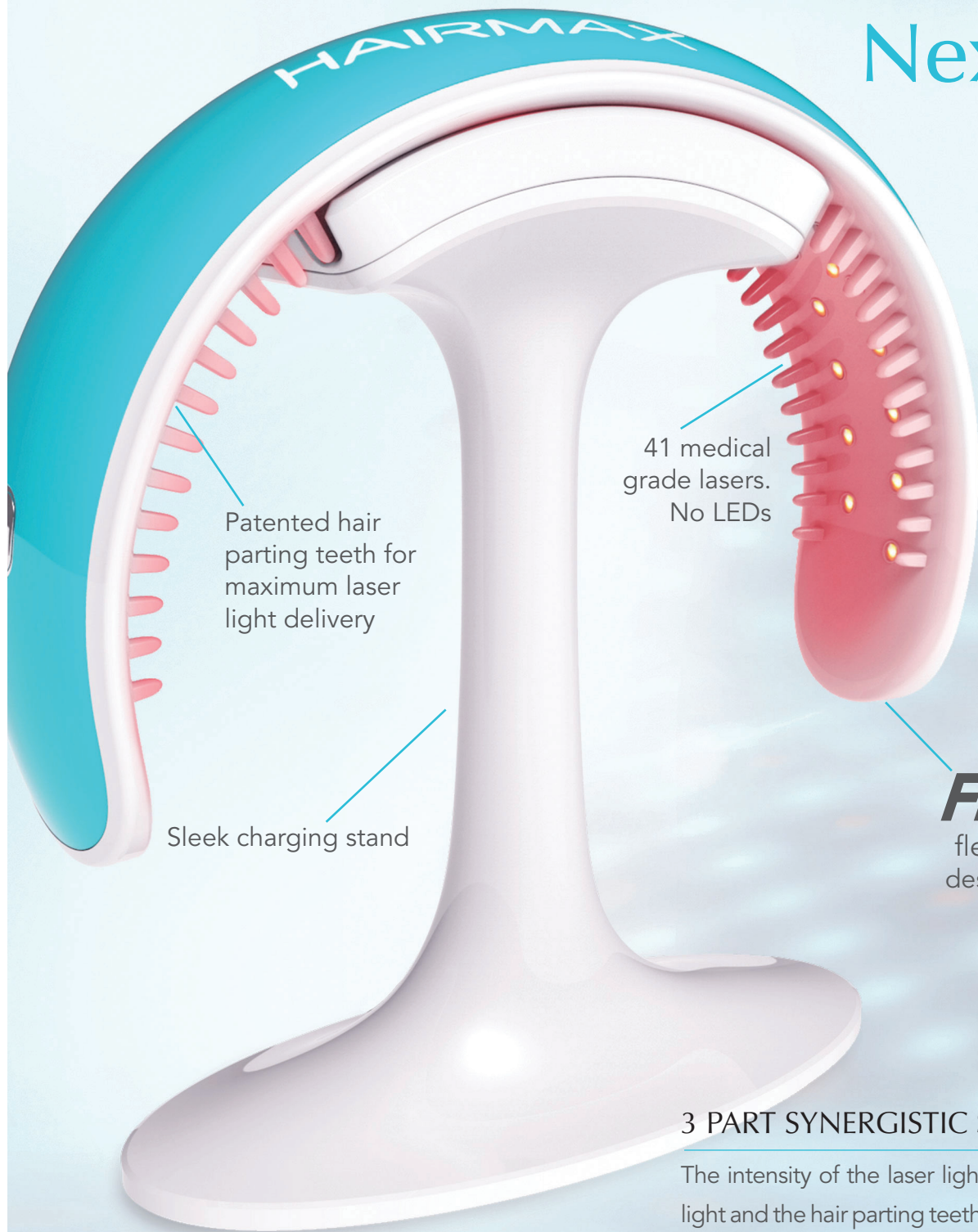

Patented hair parting teeth for maximum laser light delivery

Sleek charging stand

41 medical grade lasers. No LEDs

FLEX-FIT
flexible headband designed for comfort

CLINICALLY TESTED

FDA CLEARED*

DOCTOR RECOMMENDED

3 PART SYNERGISTIC SYSTEM

The intensity of the laser light, the density of the laser light and the hair parting teeth work in unison to provide a highly efficacious hair loss treatment.

HAIRMAX® LASERBAND 41

The LaserBand 41 is designed for maximum effectiveness, comfort and ease of use. Featuring a slim, Flex-Fit headband design that adjusts to fit all head sizes comfortably. The LaserBand delivers therapeutic light energy through 41 medical grade lasers, providing focused light energy to your hair follicles, unlike LEDs which do not. It's lightweight and battery powered, allowing for freedom of movement during treatment and portability while traveling.

SEEING IS BELIEVING



Female: Denser, fuller hair across the entire scalp.



Male: Significant new hair growth in the frontal area.

NOURISH. STIMULATE. GROW.



- Stimulates Hair Growth
- Treats Hair Loss
- Reverses the Thinning Process
- Increases Density and Fullness
- Revitalizes Damaged Hair
- No Harmful Side Effects
- Drug-Free Treatment

FDA CLEARED*
FOR MEN & WOMEN



*HairMax laser devices are indicated to treat Androgenetic Alopecia, and promote hair growth in males who have Norwood Hamilton Classifications of IIa to V and in females who have Ludwig (Savin) I-4, II-1, II-2, or frontal patterns of hair loss and who both have Fitzpatrick Skin Types I to IV. Laser Light - Avoid direct eye exposure.